

980nm Fiber Coupled DFB Laser Source

400-600mW, 2nm, Benchtop or Module



DATASHEET

[Return to the Webpage](#)



Features

- Turnkey Laser Source
- High Stability
- Advanced Feedback Control

Applications

- Medical Laser Treatment
- Biotechnology
- Others



Agiltron provides cost-effective fiber-coupled laser sources with a wide range emitting spectrum from 370nm to 2000nm and line width from 10kHz to broadband to select. Each benchtop laser source features a pigtailed laser and high-precision, low-noise auto-feedback drive electronics to ensure constant output power or a constant driving current, and an integrated temperature control unit maintains optimal operating conditions. Each unit features a front fiber output connector and a universal power supply compatible with 100 to 240 VAC. We offer two packages: benchtop for ease of use and compact module for system integration. The user interface benchtop includes an intuitive LCD display for independent control of output power and temperature via two front rotating knobs. The module has two front output power and temperature settings. All units have a built-in isolator option to prevent reflection-induced laser emissions instability. We produce fiber-coupled isolators from 370nm to 2000nm. An isolator is essential to obtain stable laser output.

Specifications

Parameter	Min	Typical	Max	Unit
Operating Case Temperature	-5	-	75	°C
Storage Temperature	-40	-	85	°C
LD Forward Current	-	-	1100	mA
LD Reverse Current	-	-	10	μA
LD Reverse Voltage	-	-	2	V
PD Forward Current	-	-	-10	mA
PD Reverse Voltage	-	-	20	V
TEC Current	-	-	2	A
TEC Voltage	-	-	3.5	V
Fiber Bend Radius	30	-	-	mm
Relative Humidity	0	-	95%	-
Lead Soldering Time	-	-	10	s
Fiber Axial Pull Force	-	-	5	N
Fiber Side Pull Force	-	-	2.5	N
Optical & Electrical Characteristics(Tc=25 °C, unless otherwise noted)				
LD Threshold Current	-	90	110	mA
Outpower	-	-	400	mW
LD Forward Current	-	-	900	mA
Kink Free Power	450	-	-	mW
Kink Free Current	≥ 1.2*If(BOL)			mA
LD Forward Voltage	-	-	2.5	V
Center Wavelength	975	976	977	nm
Peak Wavelength Turning	-	-	0.02	nm/°C
Spectrum Width	-	-	2	nm
Spectrum Stability	-0.5	-	0.5	nm
Monitor Responsivity	-	1	20	uA/mW
Monitor Responsivity Stability	-	-	20%	-
Monitor Dark Current	-	-	50	nA
TEC Current	-	-	2	A
TEC Voltage	-	-	3.5	V
TEC Modul Power Consumption	-	-	5	W
Power Stability	>20mW		0.2	dB
	10-20mW		0.5	
	3.5-10mW		1	
Tracking Error	-0.5	-	0.5	dB
Thermistor Resistance	9.5	10	10.5	Kohm
Thermistor B constant	-	3900	-	k

Note: The specifications provided are for general applications with a cost-effective approach. If you need to narrow or expand the tolerance, coverage, limit, or qualifications, please [click this link](#):

Rev 09/20/25

© Photonwares Corporation

P +1 781-935-1200

E sales@photonwares.com

W www.agiltron.com

Information contained herein is deemed to be reliable and accurate as of the issue date. Photonwares reserves the right to change the design or specifications at any time without notice. Agiltron is a registered trademark of Photonwares Corporation in the U.S. and other countries.

980nm Fiber Coupled DFB Laser Source

400-600mW, 2nm, Benchtop or Module



DATASHEET

Benchtop Laser Source Operation Manual



- Plug in power cable
- Turn on Power Switch
- Setting the Output Power by rotating the knob
- Setting the laser diode Temperature by rotating the knob
- Connect a fiber path cable with matching connector (FC/APC is the default)
- Push the Emission switch to turn on the laser
- Measure the output power to verify

Module Laser Source Operation Manual



- Plug in power cable
- Turn on Power Switch
- Setting the Output Power by rotating the screw
- Setting the laser diode Temperature by rotating the screw
- Connect a fiber path cable with matching connector (FC/APC is the default)
- Push the Emission switch to turn on the laser
- Measure the output power to verify


980nm Fiber Coupled DFB Laser Source

400-600mW, 2nm, Benchtop or Module





DATASHEET

Module (mm): Plug-Play (Power Supply Included)

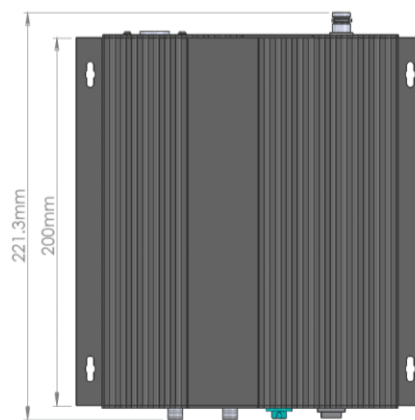

Benchtop (mm): Plug-Play 100-240VAC, USB

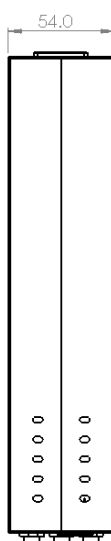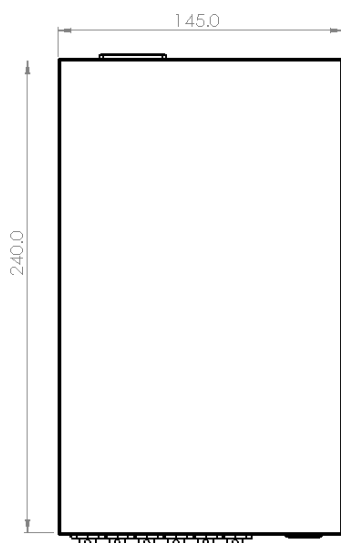

*Product dimensions may change without notice. This is sometimes required for non-standard specifications.

© Photonwares Corporation

P +1 781-935-1200

E sales@photonwares.com

W www.agiltron.com

Information contained herein is deemed to be reliable and accurate as of the issue date. Photonwares reserves the right to change the design or specifications at any time without notice. Agiltron is a registered trademark of Photonwares Corporation in the U.S. and other countries.


980nm Fiber Coupled DFB Laser Source

400-600mW, 2nm, Benchtop or Module



DATASHEET

Typical Spectrum

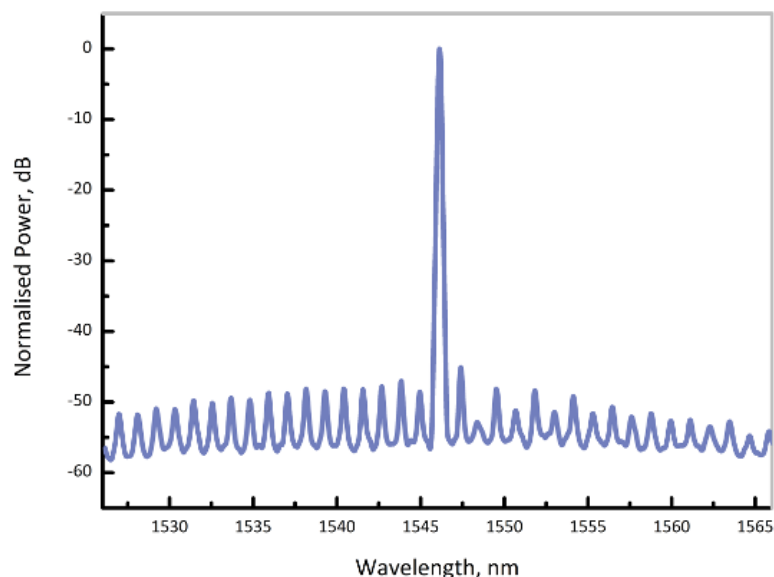

Ordering Information (Part Number)

Prefix	Wavelength	Power	Linewidth	Package	Isolator	Control Mode	TEC Cooling	Fiber Type	Connector ^[1]
FCLS-	980nm = 0980	>400mW = 4 >600mW = 6	2nm = 1	Benchtop = 1 Module = 2	None = 1 Yes = 2	Constant Current = 2 Constant Power = 1	No = 1 Yes = 2	Hi1060 = 2 PM980 = 9 50/125 = A Special = 0	FC/APC = 3 FC/PC = 2 Non = 1 SC/PC = 4 SC/APC = 5 LC/PC = 7 LC/UPC = U Special = 0

[1]. The default connector configuration uses fiber with 0.9 mm buffer protection. The connector cannot be installed directly onto bare fiber because the bare fiber is prone to damage during shipping. However, the connector can be assembled on bare fiber if a 3 cm protective loose tube is added for reinforcement. The customer may remove this protective tube after testing. The optical power handling of a standard connector is less than 0.5 W for SMF-28 fiber and decreases further for smaller-core fibers.